

medicine®

ReMed+

The smart Pillbox like never before




Contents

1. What's in the Box	2
2. Product Features	3
3. Getting to Know Your ReMed+	4
4. Quick Start	6
5. Setup ReMed+	
5.1 Configure Your Device	7
5.2 Add Medication	8
5.3 Adherence Report	10
5.4 Battery Information	11
5.5 Snooze Button	12
5.6 Reset Device	12
6. Safety Information	13
7. Warranty	14
8. Product Information	15
9. Contact Information	16

1. What's in the Box

- ReMed+
- USB-C Type Cable
- User Manual
- Thank You Card

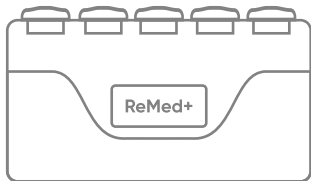

(2)

2. Product Features



Medicine
Reminders



Mobile App
Adherence
Reports



Snooze
Reminders



WiFi



USB C
Port




Light
Indication




10 Days
Battery Life

3. Getting to Know Your ReMed+

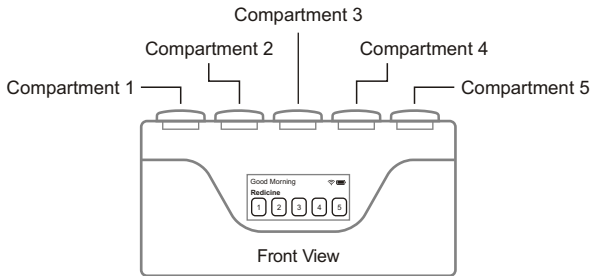

(4)

3. Getting to Know Your ReMed+

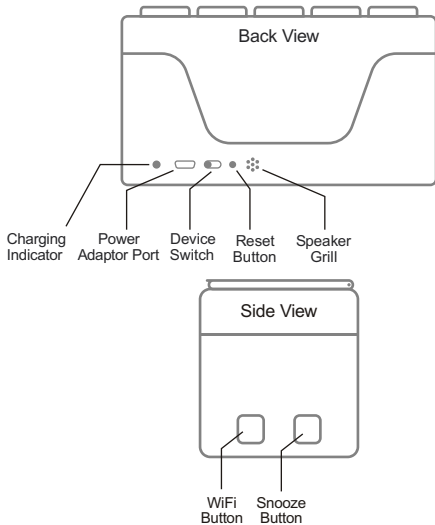

(5)

4. Quick Start

- Move the slide switch towards the USB port to turn ON the device. The display will take 5 seconds, and the home screen will be displayed.
- Tap the WiFi button. The display will show the WiFi symbol.
- Open WiFi setting on your Mobile and connect to "REMED_XXXX". To connect to the WiFi, Scan the QR code given on the sticker below the device, or enter the WiFi credentials provided on it. Now, open the Medicine Application and schedule your Medications.

5. Using ReMed+


5.1 Configure Your Device

- Connect the Device to the mobile via WiFi.
- Open Redicine Application and tap on "WiFi" symbol.
- There you will see your device in the connected device.
- Click on that connected device and you will navigate to "Pillbox Settings"
- Click on "Name" and change your device name.
(Character limit 15)
- You can see the battery percentage on this page.
- You can turn ON/OFF the device "Alarm volume" from the application.
- You can change the "Snooze Duration" time from the application. When you press the Snooze Button on the device, it will remind you to take your current medication after the time set in Snooze Duration.

5.2 Add Medication

- Connect the Device to mobile via WiFi.
- Open Redicine Application and tap on +
- Fill in all the details of your medication information.
- Using the application, select a **compartment** that is suitable for your medication, and place the medication inside the selected compartment.
- For **How often do you take?** you can select an option based on how frequently you have to take medications.
 - If you are required to consume medication on a daily basis, you can select the option of "Everyday".
 - If you are required to consume medications on some selected days of the week, tap on "Selected days".

(8)

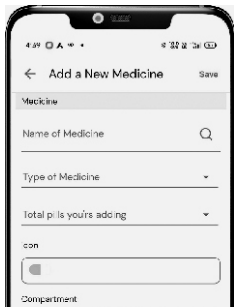

5.2 Add Medication


- If you are required to consume medication on some days interval, tap on "Days interval". For example, if you have to take the medication once a day in the thrice a day interval, you will have to select this option and 3 option from the slider.
- **How many times a day** indicates how frequently medicines consumed in a day. For example: If a doctor prescribes a medication to be taken "Three times a day", then select that option from the application and add dosage and timing according to the prescription.
- **Smart light window** lets you set the alarm with unique variations. For example: If your usual mealtime is between 8:00 pm to 8:30 pm, you can set a half-hour or hour of light window to conveniently consume your medicine. The smart alarm window allows you to consume your medicines conveniently on time, and keep a healthy record of your medication adherence.

5.3 Adherence Report

- An adherence report shows how well a patient follows a treatment plan, like if they are taking their medicine as prescribed. It helps doctors/caregivers to see if patient is sticking to the prescription.

To see the adherence report in the mobile application, Please follow the below given steps :

- Connect device to the mobile via WiFi.
- Open Medicine Application and tap on Report.
- Click on Sync now and report will be loaded on your mobile screen.

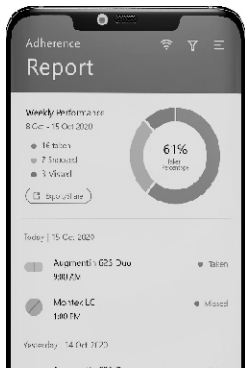

(10)

5.4 Battery Information

- The battery life of your ReMed+ is up to 10 days on a single charge.
- The battery symbol is shown on the display.
 - If the battery symbol shows empty (☐) on the display then connect the charger to the device immediately. During charging device's charging LED will be Red and charging symbol (⚡) will display.
 - When the battery symbol becomes full (■) from charging, it means that the battery is charged and the Device's charging LED will be Green.

5.5 Snooze Button

- By pressing Snooze Button once, it will snooze the medication for 10 Minutes (or according to time set from settings).
- This button also works as taking medication stored in the other compartment. To take medication from the other compartment, press and hold the Snooze Button for 5 seconds.

5.6 Reset Device

- If you ever need a fresh start with your device, the reset button is here to help. Simply press and hold it for 10 seconds. Your device will be ready for a new beginning. Whether you're passing it on to a loved one or simply want to set up it again, this feature will come in.

6. Safety Information

- Read the user manual and follow all instructions and warnings.
- Keep the product away from water and other liquids.
- Keep the product away from direct heat sources.
- Do not use the product if it appears damaged or broken.
- Do not use the product in a way that could be dangerous, such as while driving or operating heavy machinery.
- Keep the product out of reach of children and pets.
- Do not attempt to open or repair the device yourself, as this may damage the device and void the warranty.
- Notify your healthcare provider immediately if you suspect any malfunctions or security breaches.
- Use only the cable provided with the device for charging device to ensure proper functioning and safety.
- Always consult your healthcare provider before making any changes to your medication schedule.


7. Warranty

Our medicine reminder pillbox comes with a 1-year warranty for your peace of mind.

- The consumer should read the warranty information carefully and keep the receipt or proof of purchase as it may be required to make a claim.
- It's important to note that this warranty may have some limitations or exclusions, such as not covering damage caused by misuse or neglect.
- Customer service is available for any questions regarding the warranty and the product.
- The warranty will be valid for the original owner and not transferable to subsequent owners.

"To claim warranty, please contact support@redicinemedsol.com"

8. Product Information

- WiFi Network: Supports 802.11 b/g/n
- Power Input: 5V  1A (USB-C Type)
- Dimension: 215mm x 113mm x 113mm (LxWxH)
- Battery: 3.7V 1800mAh (Rechargeable)
- Display: E-Paper Display - 2.9 inch
296x128 Pixel (DPI - 125)

9. Contact Information

For any other question
Reach us

+91-9601071614

info@redicinemedsol.com

www.redicinemedsol.com

We wish you a beautiful experience with Redicine.

